**en**  
**PROCESSING (CLEANING, DISINFECTION AND STERILIZATION)**  
**OF KURZ INSTRUMENTS**

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**en**

## **PROCESSING (CLEANING, DISINFECTION AND STERILIZATION) OF KURZ INSTRUMENTS**

### **FUNDAMENTAL POINTS**

All instruments are to be cleaned, disinfected and sterilized before use (Cleaning and disinfection after removal of the protective packaging). This is required prior to each application as well as for the first use after delivery of the unsterile instruments. Effective cleaning and disinfection are an indispensable requirement for an effective sterilization.

Please ensure that only sufficiently device- and product- specific validated procedures are used for cleaning, disinfection, and sterilization, that the used devices (Washer-Disinfector (WD), sterilizer) are maintained and checked regularly and that the validated parameters are applied in each cycle.

To avoid unnecessary contamination of the sterilization-tray, please collect contaminated instruments separately (without putting them back into the sterilization tray). Pre-clean, clean, disinfect and check the instruments as well as the sterilization tray before sorting the instruments back into the sterilization tray. The completely equipped sterilization tray is then sterilized.

The sterilization trays must be cleaned and disinfected in un-loaded condition, with the lid removed and both, tray and lid facing downwards. A cleaning and disinfection of trays loaded with instruments must not be performed.

Additionally, please pay attention to the legal provisions valid for your country as well as to the hygienic instructions of the hospital or doctor's practice. This applies particularly to the different guidelines regarding the inactivation of prions (not relevant for USA), which can require the application of cleaning detergents with proven prion efficiency as well as a sterilization with more intensive parameters.

Attention: In case of some instruments additional or deviating procedures are required (see chapter „Specific Aspects“).

## INSTRUCTIONS FOR CLEANING AND DISINFECTION

### Basics

If possible, an automated procedure (WD) should be used for cleaning and disinfection of the instruments. A manual procedure – even in case of application of an ultrasonic bath – should only be used if an automated procedure is not available; in this case, the significantly lower efficiency and reproducibility of a manual procedure has to be considered.

The pre-treatment step is to be performed in both cases.

### Pre-treatment

Please remove coarse impurities of the instruments directly after application (within a maximum of 2 h).

### Procedure

1. Disassemble the instruments as far as possible (see chapter „Specific Aspects“).
2. Rinse the disassembled instruments at least 1 min under running water (temperature < 35 °C/95 °F). Agitate movable parts at least three times during pre-rinsing. If applicable (see chapter „Specific Aspects“): Rinse all lumens of the instruments at least three times with a syringe (minimum volume see chapter „Specific Aspects“).
3. Soak the instruments for the given soaking time in the pre-cleaning solution<sup>1</sup> (ultrasonic bath, ultrasound not activated) that the instruments are completely covered with the pre-cleaning solution. Pay attention that there is no contact between the instruments. Assist cleaning by careful brushing of all inner and outer surfaces with a soft brush (at beginning of soaking, aids see chapter „Specific Aspects“). Agitate movable parts at least three times during pre-cleaning.
4. If applicable (see chapter „Specific Aspects“): Rinse all lumens of the instruments at least three times at the beginning and at the end of the soaking time with a syringe (minimum volume see chapter „Specific Aspects“).
5. Activate ultrasound for an additional soaking time (but not less than 5 min).
6. Then, remove the instruments of the pre-cleaning solution and post-rinse them at least three times intensively (at least 1 min) with water. Agitate movable parts at least three times during post-rinsing.
7. If applicable (see chapter „Specific Aspects“): Rinse all lumens of the instruments at least three times at the beginning and at the end of the soaking time with a syringe (minimum volume see chapter „Specific Aspects“).

Pay attention to the following points during selection of the cleaning detergent<sup>1</sup>:

- fundamental suitability for the cleaning of instruments made of metallic or plastic material.
- suitability of the cleaning detergent for ultrasonic cleaning (no foam development).
- compatibility of the cleaning detergent with the instruments (see chapter „material resistance“).

Pay attention to the instructions of the detergent manufacturer regarding concentration, temperature and soaking time as well as post-rinsing. Please use only freshly prepared solutions as well as only sterile or low contaminated water (max. 10 germs/ml) as well as low endotoxin contaminated water (max. 0.25 endotoxin units/ml), for example purified/highly purified water, and a soft, clean, and lint-free cloth and/or filtered air for drying, respectively.

<sup>1</sup> In case of application of a cleaning and disinfection detergent for the pre-cleaning (e.g. in consequence of personnel's safety) please consider, that the pre-cleaning solvent should be aldehyde-free (otherwise blood impurities), possess a fundamentally approved efficiency (for example VAH/DGHM or FDA/EPA approval/clearance/registration or CE marking), be suitable for the disinfection of instruments made of metallic or plastic material, and be compatible with the instruments (see chapter „material resistance“). Please consider, that a disinfectant used in the pre-treatment step serves only the personnel's safety, but cannot replace the disinfection step later to be performed after cleaning.

### AUTOMATED CLEANING/DISINFECTION (WD)

Pay attention to the following points during selection of the WD:

- WD according to EN ISO/ANSI AAMI ST15883 and with fundamentally approved efficiency (for example CE marking according to EN ISO 15883 or DGHM or FDA approval/clearance/registration).
- if possible, selection of an approved WD-program for thermal disinfection (AO value ≥ 3000 or – in case of older devices – at least 5 min at 90°C/194°F; in case of chemical disinfection danger of remnants of the disinfectant on the instruments),
- fundamental suitability of the WD-program for instruments,
- WD-program with at least three degrading steps after cleaning and a sufficient number of rinsing steps (respectively neutralization, if applied) or contact-based rinsing control recommended in order to prevent effectively remnants of the detergents),
- post-rinsing only with sterile or low contaminated water (max. 10 germs/ml, max. 0.25 endotoxin units/ml), for example purified/highly purified water,
- only use of filtered air (oil-free, low contamination with microorganisms and particles) for drying,
- regularly maintenance and check/calibration of the WD.

Pay attention to the following points during selection of the cleaning detergent:

- fundamental suitability for the cleaning of instruments made of metallic or plastic material
- additional application – in case of non-application of a thermal disinfection – of a suitable disinfectant with approved efficiency (for example VAH/ DGHM or FDA/ EPA approval/clearance/registration or CE marking) compatible to the used cleaning detergent
- compatibility of the used detergents with the instruments (see chapter „material resistance“)

Pay attention to the instructions of the detergent manufacturers regarding concentration, temperature and soaking time as well as post-rinsing.

### Procedure

1. Disassemble the instruments as far as possible (see chapter „Specific Aspects“).
2. Transfer the disassembled instruments in the WD (pay attention that the instruments have no contact). If applicable (see chapter „Specific Aspects“): Connect the instruments to the rinsing port of the WD.
3. Start the program.
4. Disconnect and remove the instruments of the WD after end of the program.
5. Check and pack the instruments immediately after the removal (see chapters „check“, „maintenance“, and „packaging“, if necessary, after additional post-drying at a clean place).

The fundamental suitability of the instruments for an effective automated cleaning and disinfection was demonstrated by an independent, governmentally accredited and recognized (§ 15 (5) MPG) test laboratory by application of an ultrasonic bath of the SONOREX series with 35 kHz (BANDELIN electronic, Berlin) for pre-cleaning, the WD G 7836 CD, Miele & Cie. GmbH & Co., Gütersloh (thermal disinfection) and the pre-cleaning and cleaning detergent Neodisher medizym or Neodisher mediclean forte (both Dr. Weigert GmbH & Co. KG, Hamburg) considering the specified procedure.

## MANUAL CLEANING AND DISINFECTION

Pay attention to the following points during selection of the cleaning and disinfection detergents:

- fundamental suitability for the cleaning and disinfection of instruments made of metallic or plastic material
- suitability of the cleaning detergent for ultrasonic cleaning (no foam development)
- application of a disinfectant with approved efficiency (for example VAH/DGHM or FDA/EPA approval/clearance/registration or CE marking) compatible with the used cleaning detergent
- compatibility of the used detergents with the instruments (see chapter „material resistance“)

Combined cleaning/disinfection detergents should not be used, except in case of extremely low contamination (no visible impurities, not applicable for USA).

Pay attention to the instructions of the detergent manufacturers regarding concentration, temperature and soaking time as well as post-rinsing. Please use only freshly prepared solutions as well as only sterile or low contaminated water (max. 10 germs/ml) as well as low endotoxin contaminated water (max. 0.25 endotoxin units/ml), for example purified/highly purified water, and a soft, clean, and lint-free cloth and/or filtered air for drying, respectively.

### Procedure

#### Cleaning

1. Disassemble the instruments as far as possible (see chapter „Specific Aspects“).
2. Soak the disassembled instruments for the given soaking time in the cleaning solution (ultrasonic bath, ultrasound not activated so that the instruments are completely covered. Pay attention that there is no contact between the instruments. Assist cleaning by careful brushing of all inner and outer surfaces (at beginning of soaking, aids see chapter „Specific Aspects“). Agitate movable parts at least three times during cleaning.
3. If applicable (see chapter „Specific Aspects“): Rinse all lumens of the instruments at least five times at the beginning and at the end of the soaking time by application of a single-use syringe (minimum volume see chapter „Specific Aspects“).
4. Activate ultrasound for an additional soaking time (but not less than 5 min).
5. Then, remove the instruments of the cleaning solution and post-rinse them at least three times intensively (at least 1 min) with water. Agitate movable parts at least three times during post-rinsing.
6. If applicable (see chapter „Specific Aspects“): Rinse all lumens of the instruments at least five times at the beginning and at the end of the soaking time by application of a single-use syringe (minimum volume see chapter „Specific Aspects“).
7. Check the instruments (see chapters „check“ and „maintenance“).

### Disinfection

1. Soak the disassembled, cleaned, and checked instruments for the given soaking time in the disinfectant solution so that the instruments are completely covered. Pay attention that there is no contact between the instruments. Agitate movable parts at least three times during disinfection.
2. If applicable (see chapter „Specific Aspects“): Rinse all lumens of the instruments at least five times at the beginning and at the end of the soaking time by application of a single-use syringe (minimum volume see chapter „Specific Aspects“).
3. Then, remove the instruments of the disinfectant solution and post-rinse them at least five times intensively (at least 1 min) with water. Agitate movable parts at least three times during post-rinsing.
4. If applicable (see chapter „Specific Aspects“): Rinse all lumens of the instruments at least five times at the beginning and at the end of the soaking time by application of a single-use syringe (minimum volume see chapter „Specific Aspects“).
5. Dry by blowing off/through with filtered compressed air.
6. Pack the instruments immediately after the removal (see chapter „packaging“, if necessary after additional post-drying at a clean place).

The fundamental suitability of the instruments for an effective manual cleaning and disinfection was demonstrated by an independent, governmentally accredited and recognized (§ 15 (5) MPG) test laboratory by application of an ultrasonic bath of the SONOREX series with 35 kHz (BANDELIN electronic, Berlin) for pre-cleaning and manual cleaning as well as the pre-cleaning and cleaning detergent Cidezyme/Enzol and the disinfectant Cidex OPA (ASP Deutschland GmbH, Aldingen) considering the specified procedure and the records are on file at Heinz Kurz GmbH.

### CHECK

Check all instruments after cleaning or cleaning/disinfection, respectively, on corrosion, damaged surfaces, still readable marking, and impurities as well as free movement. Do not further use damaged instruments (for limitation of the numbers of re-use cycles see chapter „reusability“). Still dirty instruments are to be cleaned and disinfected again.

### MAINTENANCE

Assemble disassembled instruments again as displayed in chapter „Specific Aspects“. Instrument oils or grease must not be applied.

### PACKAGING

Please insert the cleaned and disinfected instruments in the corresponding sterilization trays and pack them in single-use sterilization packagings (single or double packaging), which fulfill the following requirements (material/process):

- EN ISO/ANSI AAMI ISO 11607 (for USA: FDA clearance)
- suitable for steam sterilization (temperature resistance up to at least 138°C (280°F), sufficient steam permeability)
- sufficient protection of the instruments as well as of the sterilization packagings to mechanical damage

## STERILIZATION

Please use for sterilization only the listed sterilization procedures; other sterilization procedures must not be applied.

### Steam sterilization

- fractionated vacuum/dynamic air removal procedure<sup>2,3</sup> (with sufficient product drying<sup>4</sup>)
- steam sterilizer according to EN 13060/EN 285 or ANSI AAMI ST79 (for USA: FDA clearance)
- validated according to EN ISO 17665 (valid IQ/OQ (commissioning) and product specific performance qualification (PQ))
- maximum sterilization temperature 134°C (273°F; plus tolerance according to EN ISO 17665)
- sterilization time (exposure time at the sterilization temperature):

| Area            | fractionated vacuum/dynamic air removal                          | gravity displacement |
|-----------------|------------------------------------------------------------------|----------------------|
| Germany         | 5 min at 134°C (273°F)                                           | not recommended      |
| Switzerland     | 18 min at 134°C (273°F)                                          | not recommended      |
| USA             | 4 min at 132°C (270°F), drying time at least 20 min <sup>4</sup> | not recommended      |
| other countries | at least 3 min <sup>5</sup> at 132°C (270°F) / 134°C (273°F)     | not recommended      |

2 at least three vacuum steps

3 The less effective gravity displacement procedure must not be used in case of availability of the fractionated vacuum procedure, requires significantly longer sterilization times as well as a sterilizer, procedure, parameter, and product specific process development and validation under sole responsibility of the user.

4 The effectively required drying time depends directly on parameters in sole responsibility of the user (load configuration and density, sterilizer conditions, ...) and by this is to be determined by the user. Nevertheless, drying times less than 20 min must not be applied.

5 respectively 18 min (inactivation of prions, not relevant for USA)

## STORAGE

Please store the instruments after sterilization in the sterilization packagings at a dry and dust-free place.

## MATERIAL RESISTANCE

Please take care that the listed substances are not ingredients of the cleaning or disinfection detergent:

- organic, mineral, and oxidizing acids (minimum admitted pH-value 5.5)
- stronger lyes (maximum admitted pH-value 11, neutral/enzymatic or weak alkaline cleaning detergent recommended)
- organic solvents (for example: acetone, ether, alcohol, benzine)
- oxidizing agents (for example: peroxide)
- halogens (chlorine, iodine, bromine)
- aromatic, halogenated hydrocarbons
- oils

Please consider during selection of the detergents in addition, that corrosion inhibitors, neutralizing agents, and/or rinse aids may cause potential critical remnants on the instruments. Acid neutralizing agents or rinse aids must not be applied.

Please do not clean any instruments and sterilization trays by use of metal brushes or steel wool.

Please do not expose any instruments and sterilization trays to temperatures higher than 138°C (280°F)!

## REUSABILITY

The instruments can be reused – in case of adequate care and if they are undamaged and clean. The user is responsible for each further use as well as for the use of damaged and dirty instruments (no liability in case of disregard).

## DISPOSAL

Please refer to country-specific disposal guidelines.

|                                                                                                                |                                                                                                                                                                                                                                                                         |
|----------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>Art. no</b>                                                                                                 | 8000100 and 8000106                                                                                                                                                                                                                                                     |
| <b>article specification</b>                                                                                   | KURZ Meter<br>KURZ Meter Set (see Figure 1)                                                                                                                                                                                                                             |
| <b>rinsing volume</b>                                                                                          | Holder inside: 100 ml (single-use syringe)<br>tube inside: single-use syringe (0.1-1 ml) and needle                                                                                                                                                                     |
| <b>brush</b>                                                                                                   | Holder inside: pipe cleaner 6 mm<br>outside: standard brush/tooth brush                                                                                                                                                                                                 |
| <b>specific/additional procedure in case of</b>                                                                |                                                                                                                                                                                                                                                                         |
| pretreatment                                                                                                   | Disassemble instrument into three units: „Holder“, „Handle/pusher/probe“, „Tube“<br><br>Brush outside and inside.<br>Flush blind bore and tube<br>Check the blind bore for leftovers.                                                                                   |
| manual cleaning/ disinfection                                                                                  | Disassembled.<br><br>Outside: standard brush, holder inside: pipe cleaner (6 mm) and rinse syringe single-use (50–100 ml), tube inside: single-use syringe (0.1-1 ml) and needle, US treatment.<br><br>Do not apply specific instrument tray for cleaning/disinfection! |
| automated cleaning/ disinfection                                                                               | Disassembled.<br>Small pieces basket. Do not apply specific instrument tray for cleaning/disinfection                                                                                                                                                                   |
| maintenance/ assembly                                                                                          | No assembly prior to sterilization<br>Lubrication not admitted!                                                                                                                                                                                                         |
| <b>packaging</b>                                                                                               | Disassembled in specific instrument tray, paper/film bag (single or twice), lubrication prohibited                                                                                                                                                                      |
| <b>sterilization</b>                                                                                           | Disassembled in specific sterilization tray.                                                                                                                                                                                                                            |
| recommended classification according to KRINKO/RKI/BfArM guidance (only Germany, with respect to intended use) | Critical B                                                                                                                                                                                                                                                              |

|                                                                                                                |                                                                                                                                                                         |
|----------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>Art. no</b>                                                                                                 | 8000108                                                                                                                                                                 |
| <b>article specification</b>                                                                                   | Ventilation Tube Insertion Forceps                                                                                                                                      |
| <b>rinsing volume</b>                                                                                          | -                                                                                                                                                                       |
| <b>brush</b>                                                                                                   | standard brush/tooth brush                                                                                                                                              |
| <b>specific/additional procedure in case of</b>                                                                |                                                                                                                                                                         |
| pretreatment                                                                                                   | Open the latch on top of the pliers and open the pliers.<br>Brush. Agitate moving parts during soaking and rinsing.<br>Open pliers to maximum for ultrasonic treatment. |
| manual cleaning/ disinfection                                                                                  | With opened latch.<br><br>Brush.<br><br>Agitate moving parts during soaking and rinsing.<br>Open pliers to maximum for ultrasonic treatment.                            |
| automated cleaning/ disinfection                                                                               | Small pieces basket                                                                                                                                                     |
| maintenance/ assembly                                                                                          | Lubrication not admitted!                                                                                                                                               |
| <b>packaging</b>                                                                                               | Paper/film bag (single or twice).                                                                                                                                       |
| <b>sterilization</b>                                                                                           | Close the latch. Not lubricated.                                                                                                                                        |
| recommended classification according to KRINKO/RKI/BfArM guidance (only Germany, with respect to intended use) | Critical B                                                                                                                                                              |

|                                                                                                                                 |                                                                                                  |
|---------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------|
| <b>Art. no</b>                                                                                                                  | 8000109                                                                                          |
| <b>article specification</b>                                                                                                    | Malleus Handle Cavity Bending Pliers                                                             |
| <b>rinsing volume</b>                                                                                                           | -                                                                                                |
| <b>brush</b>                                                                                                                    | standard brush/tooth brush                                                                       |
| <b>specific/additional procedure in case of</b>                                                                                 |                                                                                                  |
| pretreatment                                                                                                                    | Brush.<br>Agitate moving parts during soaking and rinsing.<br>Check pointed gap for remnants.    |
| manual cleaning/<br>disinfection                                                                                                | Brush.<br>Agitate moving parts<br>during soaking and rinsing. Check pointed gap<br>for remnants. |
| automated<br>cleaning/<br>disinfection                                                                                          | Small pieces basket                                                                              |
| maintenance/<br>assembly                                                                                                        | Lubrication not admitted!                                                                        |
| <b>packaging</b>                                                                                                                | Paper/film bag (single or twice).                                                                |
| <b>sterilization</b>                                                                                                            | Not lubricated!                                                                                  |
| recommended<br>classification according<br>to KRINKO/RKI/BfArM<br>guidance (only Germa-<br>ny, with respect to<br>intended use) | Critical B                                                                                       |

|                                                                                                                                 |                                   |
|---------------------------------------------------------------------------------------------------------------------------------|-----------------------------------|
| <b>Art. no</b>                                                                                                                  | 8000127                           |
| <b>article specification</b>                                                                                                    | Soft-Clip Hook                    |
| <b>rinsing volume</b>                                                                                                           | -                                 |
| <b>brush</b>                                                                                                                    | standard brush/tooth brush        |
| <b>specific/additional procedure in case of</b>                                                                                 |                                   |
| pretreatment                                                                                                                    | Brush.                            |
| manual cleaning/<br>disinfection                                                                                                | Brush.                            |
| automated<br>cleaning/<br>disinfection                                                                                          | Small pieces basket               |
| maintenance/<br>assembly                                                                                                        | Lubrication not admitted!         |
| <b>packaging</b>                                                                                                                | Paper/film bag (single or twice). |
| <b>sterilization</b>                                                                                                            | Not lubricated!                   |
| recommended<br>classification according<br>to KRINKO/RKI/BfArM<br>guidance (only Germa-<br>ny, with respect to<br>intended use) | Critical B                        |

|                                                                                                                                 |                                                                                                                |
|---------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------|
| <b>Art. no</b>                                                                                                                  | 8000131                                                                                                        |
| <b>article specification</b>                                                                                                    | TVFMI<br>Bending Pliers                                                                                        |
| <b>rinsing volume</b>                                                                                                           | -                                                                                                              |
| <b>brush</b>                                                                                                                    | standard brush/tooth brush                                                                                     |
| <b>specific/additional procedure in case of</b>                                                                                 |                                                                                                                |
| pretreatment                                                                                                                    | Brush.<br>Agitate moving parts during soaking and rinsing.<br>Open pliers to maximum for ultrasonic treatment. |
| manual cleaning/<br>disinfection                                                                                                | Brush.<br>Agitate moving parts during soaking and rinsing.<br>Open pliers to maximum for ultrasonic treatment. |
| automated<br>cleaning/<br>disinfection                                                                                          | Small pieces basket                                                                                            |
| maintenance/<br>assembly                                                                                                        | Lubrication not admitted!                                                                                      |
| <b>packaging</b>                                                                                                                | Paper/film bag (single or twice).                                                                              |
| <b>sterilization</b>                                                                                                            | Closed. Not lubricated!                                                                                        |
| recommended<br>classification according<br>to KRINKO/RKI/BfArM<br>guidance (only Germa-<br>ny, with respect to<br>intended use) | Critical B                                                                                                     |

|                                                                                                                                 |                                                                                                                                                                                                                        |
|---------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>Art. no</b>                                                                                                                  | 8000135                                                                                                                                                                                                                |
| <b>article specification</b>                                                                                                    | Vario Holder (see figure 2)                                                                                                                                                                                            |
| <b>rinsing volume</b>                                                                                                           | 50-100 ml (single-use syringe)                                                                                                                                                                                         |
| <b>brush</b>                                                                                                                    | Vario-Holder inside: pipe cleaner 6 mm<br>outside: standard brush/tooth brush)                                                                                                                                         |
| <b>specific/additional procedure in case of</b>                                                                                 |                                                                                                                                                                                                                        |
| pretreatment                                                                                                                    | Disassemble into three units: Upper Part, Lower Part, Pin<br>Outside: standard brush, Inner diameter: pipe cleaner (6 mm) and rinse single-use syringe (50-100 ml)<br><br>Do not insert into specific instrument tray! |
| manual cleaning/<br>disinfection                                                                                                | Outside: standard brush, Inner diameter: pipe cleaner (6 mm) and rinse single-use syringe (50-100 ml)<br><br>Brush.<br><br>Do not apply specific instrument tray for cleaning/disinfection!                            |
| automated<br>cleaning/<br>disinfection                                                                                          | Small pieces basket, Check pointed gap for remnants.<br><br>Do not apply specific instrument tray for cleaning/disinfection!                                                                                           |
| maintenance/<br>assembly                                                                                                        | No assembly prior to sterilization<br>Lubrication not admitted!                                                                                                                                                        |
| <b>packaging</b>                                                                                                                | Disassembled in specific instrument tray, paper/film bag (single or twice), lubrication prohibited.                                                                                                                    |
| <b>sterilization</b>                                                                                                            | Disassembled in specific sterilization tray.                                                                                                                                                                           |
| recommended<br>classification according<br>to KRINKO/RKI/BfArM<br>guidance (only Germa-<br>ny, with respect to<br>intended use) | Uncritical                                                                                                                                                                                                             |

|                                                                                                                                 |                                                                                                 |
|---------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------|
| <b>Art. no</b>                                                                                                                  | 8000136                                                                                         |
| <b>article specification</b>                                                                                                    | Titanium Tweezers                                                                               |
| <b>rinsing volume</b>                                                                                                           | -                                                                                               |
| <b>brush</b>                                                                                                                    | standard brush/tooth brush                                                                      |
| <b>specific/additional procedure in case of</b>                                                                                 |                                                                                                 |
| pretreatment                                                                                                                    | Brush.<br>Check pointed gap for remnants.                                                       |
| manual cleaning/<br>disinfection                                                                                                | Brush.<br>Check pointed gap for remnants.                                                       |
| automated<br>cleaning/<br>disinfection                                                                                          | Small pieces basket,<br><br>Do not apply specific instrument tray<br>for cleaning/disinfection! |
| maintenance/<br>assembly                                                                                                        | Lubrication not admitted                                                                        |
| <b>packaging</b>                                                                                                                | Paper/film bag (single or twice).                                                               |
| <b>sterilization</b>                                                                                                            | Closed.<br>Not lubricated                                                                       |
| recommended<br>classification according<br>to KRINKO/RKI/BfArM<br>guidance (only Germa-<br>ny, with respect to<br>intended use) | Uncritical                                                                                      |

|                                                                                                                                 |                                                                                                                                                 |
|---------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>Art. no</b>                                                                                                                  | 8000137                                                                                                                                         |
| <b>article specification</b>                                                                                                    | Titanium Micro Closing Forceps                                                                                                                  |
| <b>rinsing volume</b>                                                                                                           | -                                                                                                                                               |
| <b>brush</b>                                                                                                                    | standard brush/tooth brush                                                                                                                      |
| <b>specific/additional procedure in case of</b>                                                                                 |                                                                                                                                                 |
| pretreatment                                                                                                                    | Unhook spring and open up to 90°.<br>Brush.<br><br>Agitate moving parts during soaking and rinsing.<br>Open up to 90° for ultrasonic treatment. |
| manual cleaning/<br>disinfection                                                                                                | Spring unhooked.<br>Brush.<br>Agitate moving parts during soaking and rinsing.<br>Open up to 90° for ultrasonic treatment.                      |
| automated<br>cleaning/<br>disinfection                                                                                          | Small pieces basket,<br><br>Do not apply specific instrument tray for<br>cleaning/disinfection!                                                 |
| maintenance/<br>assembly                                                                                                        | Hook spring.<br>Lubrication not admitted.                                                                                                       |
| <b>packaging</b>                                                                                                                | Paper/film bag (single or twice).                                                                                                               |
| <b>sterilization</b>                                                                                                            | Closed and with hooked spring. Not lubricated                                                                                                   |
| recommended<br>classification according<br>to KRINKO/RKI/BfArM<br>guidance (only Germa-<br>ny, with respect to<br>intended use) | Uncritical                                                                                                                                      |

|                                                                                                                                 |                                                                   |
|---------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------|
| <b>Art. no</b>                                                                                                                  | 8000143                                                           |
| <b>article specification</b>                                                                                                    | Trocar Handle                                                     |
| <b>rinsing volume</b>                                                                                                           | 10 ml (single-use syringe with needle of < 0.8 mm diameter)       |
| <b>brush</b>                                                                                                                    | Brush.<br>Flush blind bore.<br>Check the blind bore for leftovers |
| <b>specific/additional procedure in case of</b>                                                                                 |                                                                   |
| pretreatment                                                                                                                    | Brush.<br>Flush blind bore.<br>Check the blind bore for leftovers |
| manual cleaning/<br>disinfection                                                                                                | Brush.<br>Flush blind bore.                                       |
| automated<br>cleaning/<br>disinfection                                                                                          | Small pieces basket                                               |
| maintenance/<br>assembly                                                                                                        | Lubrication not admitted                                          |
| <b>packaging</b>                                                                                                                | Paper/film bag (single or twice).                                 |
| <b>sterilization</b>                                                                                                            | Not lubricated.                                                   |
| recommended<br>classification according<br>to KRINKO/RKI/BfArM<br>guidance (only Germa-<br>ny, with respect to<br>intended use) | Critical B                                                        |

|                                                                                                                                 |                                                                                                                                                                                                      |
|---------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>Art. no</b>                                                                                                                  | 8000177/ 8000145/ 8000124                                                                                                                                                                            |
| <b>article specification</b>                                                                                                    | Tray KURZ Precise incl. Plate POM +<br>Tray stainless steel                                                                                                                                          |
| <b>rinsing volume</b>                                                                                                           | -                                                                                                                                                                                                    |
| <b>brush</b>                                                                                                                    | standard brush/tooth brush                                                                                                                                                                           |
| <b>specific/additional procedure in case of</b>                                                                                 |                                                                                                                                                                                                      |
| pretreatment                                                                                                                    | Remove all<br>instruments.<br><br>Dismantled<br><br>Brush                                                                                                                                            |
| manual cleaning/<br>disinfection                                                                                                | Disassembled<br><br>Brush.                                                                                                                                                                           |
| automated<br>cleaning/<br>disinfection                                                                                          | In empty state.Disassembled.<br>Standard Tray: Open side of basket and lid downwards<br><br>POM Plate: Open side downwards<br><br>Stainless Steel Box: Open side downwards in<br>small pieces basket |
| maintenance/<br>assembly                                                                                                        | Lubrication not admitted                                                                                                                                                                             |
| <b>packaging</b>                                                                                                                | Paper/film bag (single or twice)                                                                                                                                                                     |
| <b>sterilization</b>                                                                                                            | Closed. Not lubricated                                                                                                                                                                               |
| recommended<br>classification according<br>to KRINKO/RKI/BfArM<br>guidance (only Germa-<br>ny, with respect to<br>intended use) | Critical B                                                                                                                                                                                           |

|                                                 |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     |
|-------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>Art. no</b>                                  | 8000155                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                             |
| <b>article specification</b>                    | KURZ Precise Cartilage Knife Set (see Figure 3)                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     |
| <b>rinsing volume</b>                           | 50-100 ml (single-use syringe)                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      |
| <b>brush</b>                                    | standard brush/tooth brush                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                          |
| <b>specific/additional procedure in case of</b> |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     |
| pretreatment                                    | Blade holder outside and inside: standard brush/tooth brush)<br>Cutting Block inside: pipe cleaner 9-10 mm<br>Nut of Cutting Block/Cutting Block inside: pipe cleaner 7-8 mm<br>outside: standard brush/tooth brush)<br>Distance plates outside and inside: standard brush/tooth brush)                                                                                                                                                                                                                                                                                                             |
| manual cleaning/disinfection                    | Blade holder: Dismantle into lower part, upper part, and pin.<br>1. Undo the Screw.<br>2. Disassemble the blade holder and the cutting block.<br>3. Remove and dispose single-use blade (no repeated use admitted).<br>Brush outside and inside. Flush inside.<br><br>Cutting Block: Dismantle into lower and upper part, nut, and spacer plates (see picture in annex). Brush inside and outside. Flush inside.<br><br>Distance plates: Brush.<br>Check for remaining bristles (otherwise repeat pre-cleaning or dispose).<br><br>Do not apply specific instrument tray for cleaning/disinfection! |
| automated cleaning/disinfection                 | Disassembled. Lower, part, upper part and screw of Blade holder: Brush inside and outside.<br>Flush inside.<br>Do not apply specific instrument tray for cleaning/disinfection!                                                                                                                                                                                                                                                                                                                                                                                                                     |
| maintenance/assembly                            | Insert pin of Blade holder into the hole (but not into the thread!).<br>Insert the distance plates into the Stainless Steel Box. No further assembly prior to sterilization.<br>Lubrication not admitted!                                                                                                                                                                                                                                                                                                                                                                                           |

|                                                                                                                  |                                                                                                                                                                                                                                                                                                                                                     |
|------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>packaging</b>                                                                                                 | Disassembled in specific instrument tray (see picture in the annex).<br>Specific instrument tray in paper/film bag (single or twice).                                                                                                                                                                                                               |
| <b>sterilization</b>                                                                                             | Disassembled (exceptions see Assembly) in specific instrument tray. Not lubricated. After sterilization/directly prior to application:<br>1. Screw upper and lower part of the Cutting Block by means of the nut. Insert Distance plates as required.<br>2. After inserting a sterile blade, fix the blade holder in place by tightening the screw. |
| irrecommended classification according to KRINKO/RKI/BfArM guidance (only Germany, with respect to intended use) | Critical B                                                                                                                                                                                                                                                                                                                                          |

|                                                                                                                |                                                                                                                                    |
|----------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------|
| <b>Art. no</b>                                                                                                 | 8000171                                                                                                                            |
| <b>article specification</b>                                                                                   | Cutting Forceps                                                                                                                    |
| <b>rinsing volume</b>                                                                                          | -                                                                                                                                  |
| <b>brush</b>                                                                                                   | standard brush/tooth brush                                                                                                         |
| <b>specific/additional procedure in case of</b>                                                                |                                                                                                                                    |
| pretreatment                                                                                                   | Unhook spring and open up to 90°. Brush. Agitate moving parts during soaking and rinsing. Open up to 90° for ultrasonic treatment. |
| manual cleaning/ disinfection                                                                                  | Spring unhooked. Brush. Agitate moving parts during soaking and rinsing. Open up to 90° for ultrasonic treatment.                  |
| automated cleaning/ disinfection                                                                               | Small pieces basket                                                                                                                |
| maintenance/ assembly                                                                                          | Hook spring. Lubrication not admitted.                                                                                             |
| <b>packaging</b>                                                                                               | Paper/film bag (single or twice).                                                                                                  |
| <b>sterilization</b>                                                                                           | Closed and with hooked spring. Not lubricated                                                                                      |
| recommended classification according to KRINKO/RKI/BfArM guidance (only Germany, with respect to intended use) | Uncritical                                                                                                                         |

|                                                                                                                |                                                                                                                  |
|----------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------|
| <b>Art. no</b>                                                                                                 | 8000172                                                                                                          |
| <b>article specification</b>                                                                                   | Micro Scissors                                                                                                   |
| <b>rinsing volume</b>                                                                                          | -                                                                                                                |
| <b>brush</b>                                                                                                   | standard brush/tooth brush                                                                                       |
| <b>specific/additional procedure in case of</b>                                                                |                                                                                                                  |
| pretreatment                                                                                                   | Open up to 90°. Brush. Agitate moving parts during soaking and rinsing. Open up to 90° for ultrasonic treatment. |
| manual cleaning/ disinfection                                                                                  | Opened. Brush. Agitate moving parts during soaking and rinsing. Open up to 90° for ultrasonic treatment.         |
| automated cleaning/ disinfection                                                                               | Small pieces basket                                                                                              |
| maintenance/ assembly                                                                                          | Lubrication not admitted.                                                                                        |
| <b>packaging</b>                                                                                               | Paper/film bag (single or twice).                                                                                |
| <b>sterilization</b>                                                                                           | Closed and with hooked spring. Not lubricated                                                                    |
| recommended classification according to KRINKO/RKI/BfArM guidance (only Germany, with respect to intended use) | Uncritical                                                                                                       |

|                                                                                                                                 |                                                  |
|---------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------|
| <b>Art. no</b>                                                                                                                  | 8000173                                          |
| <b>article specification</b>                                                                                                    | Tray TTP-VARIAC                                  |
| <b>rinsing volume</b>                                                                                                           | -                                                |
| <b>brush</b>                                                                                                                    | standard brush/tooth brush                       |
| <b>specific/additional procedure in case of</b>                                                                                 |                                                  |
| pretreatment                                                                                                                    | Remove all instruments<br>Brush.                 |
| manual cleaning/<br>disinfection                                                                                                | Brush.                                           |
| automated<br>cleaning/<br>disinfection                                                                                          | In empty state,<br>Open side of basket downwards |
| maintenance/<br>assembly                                                                                                        | Lubrication not admitted.                        |
| <b>packaging</b>                                                                                                                | Paper/film bag (single or twice).                |
| <b>sterilization</b>                                                                                                            | Closed.<br>Not lubricated.                       |
| recommended<br>classification according<br>to KRINKO/RKI/BfArM<br>guidance (only Germa-<br>ny, with respect to<br>intended use) | Critical B                                       |

|                                                                                                                                 |                                                  |
|---------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------|
| <b>Art. no</b>                                                                                                                  | 8000174                                          |
| <b>article specification</b>                                                                                                    | Tray<br>KURZ-Meter                               |
| <b>rinsing volume</b>                                                                                                           | -                                                |
| <b>brush</b>                                                                                                                    | standard brush/tooth brush                       |
| <b>specific/additional procedure in case of</b>                                                                                 |                                                  |
| pretreatment                                                                                                                    | Remove all instruments<br>Brush.                 |
| manual cleaning/<br>disinfection                                                                                                | Brush.                                           |
| automated<br>cleaning/<br>disinfection                                                                                          | In empty state,<br>Open side of basket downwards |
| maintenance/<br>assembly                                                                                                        | Lubrication not admitted.                        |
| <b>packaging</b>                                                                                                                | Paper/film bag (single or twice).                |
| <b>sterilization</b>                                                                                                            | Closed.<br>Not lubricated.                       |
| recommended<br>classification according<br>to KRINKO/RKI/BfArM<br>guidance (only Germa-<br>ny, with respect to<br>intended use) | Critical B                                       |

|                                                                                                                                 |                                                  |
|---------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------|
| <b>Art. no</b>                                                                                                                  | 8000176                                          |
| <b>article specification</b>                                                                                                    | Tray Cartilage Punch Set                         |
| <b>rinsing volume</b>                                                                                                           | -                                                |
| <b>brush</b>                                                                                                                    | standard brush/tooth brush                       |
| <b>specific/additional procedure in case of</b>                                                                                 |                                                  |
| pretreatment                                                                                                                    | Remove all instruments<br>Brush.                 |
| manual cleaning/<br>disinfection                                                                                                | Brush.                                           |
| automated<br>cleaning/<br>disinfection                                                                                          | In empty state,<br>Open side of basket downwards |
| maintenance/<br>assembly                                                                                                        | Lubrication not admitted.                        |
| <b>packaging</b>                                                                                                                | Paper/film bag (single or twice).                |
| <b>sterilization</b>                                                                                                            | Closed.<br>Not lubricated.                       |
| recommended<br>classification according<br>to KRINKO/RKI/BfArM<br>guidance (only Germa-<br>ny, with respect to<br>intended use) | Critical B                                       |

|                                                                                                                                 |                                                                   |
|---------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------|
| <b>Art. no</b>                                                                                                                  | 8000177                                                           |
| <b>article specification</b>                                                                                                    | Tray for KURZ Precise                                             |
| <b>rinsing volume</b>                                                                                                           | -                                                                 |
| <b>brush</b>                                                                                                                    | standard brush/tooth brushh                                       |
| <b>specific/additional procedure in case of</b>                                                                                 |                                                                   |
| pretreatment                                                                                                                    | Remove all instruments<br>Brush.                                  |
| manual cleaning/<br>disinfection                                                                                                | In empty state,<br>Brush.                                         |
| automated<br>cleaning/<br>disinfection                                                                                          | In empty state,<br>Disassembled.<br>Open side of basket downwards |
| maintenance/<br>assembly                                                                                                        | Lubrication not admitted.                                         |
| <b>packaging</b>                                                                                                                | Paper/film bag (single or twice).                                 |
| <b>sterilization</b>                                                                                                            | Closed.<br>Not lubricated.                                        |
| recommended<br>classification according<br>to KRINKO/RKI/BfArM<br>guidance (only Germa-<br>ny, with respect to<br>intended use) | Critical B                                                        |

|                                                                                                                                 |                                                                                                                                                                            |
|---------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>Art. no</b>                                                                                                                  | 8000188                                                                                                                                                                    |
| <b>article specification</b>                                                                                                    | SteadyCrimP Forceps                                                                                                                                                        |
| <b>rinsing volume</b>                                                                                                           | -                                                                                                                                                                          |
| <b>brush</b>                                                                                                                    | standard brush/tooth brushh                                                                                                                                                |
| <b>specific/additional procedure in case of</b>                                                                                 |                                                                                                                                                                            |
| pretreatment                                                                                                                    | Open the latch on top of the pliers and open the pliers.<br>Brush.<br>Agitate moving parts during soaking and rinsing.<br>Open pliers to maximum for ultrasonic treatment. |
| manual cleaning/<br>disinfection                                                                                                | With opened latch.<br>Brush.<br>Agitate moving parts during soaking and rinsing.<br>Open pliers to maximum for ultrasonic treatment.                                       |
| automated<br>cleaning/<br>disinfection                                                                                          | Small pieces baske                                                                                                                                                         |
| maintenance/<br>assembly                                                                                                        | Lubrication not admitted.                                                                                                                                                  |
| <b>packaging</b>                                                                                                                | Paper/film bag (single or twice).                                                                                                                                          |
| <b>sterilization</b>                                                                                                            | Closed.<br>Not lubricated.                                                                                                                                                 |
| recommended<br>classification according<br>to KRINKO/RKI/BfArM<br>guidance (only Germa-<br>ny, with respect to<br>intended use) | Critical B                                                                                                                                                                 |

|                                                                                                                                 |                                                                                             |
|---------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------|
| <b>Art. no</b>                                                                                                                  | 8000193                                                                                     |
| <b>article specification</b>                                                                                                    | Cartilage Forceps Schimanski Design                                                         |
| <b>rinsing volume</b>                                                                                                           | -                                                                                           |
| <b>brush</b>                                                                                                                    | standard brush/tooth brushh                                                                 |
| <b>specific/additional procedure in case of</b>                                                                                 |                                                                                             |
| pretreatment                                                                                                                    | Brush.<br>Open and close during soaking and rinsing.<br><br>Check pointed gap for remnants. |
| manual cleaning/<br>disinfection                                                                                                | Brush.<br>Open and close during soaking and rinsing.<br><br>Check pointed gap for remnants. |
| automated<br>cleaning/<br>disinfection                                                                                          | Small pieces basket.<br><br>Check pointed gap for remnants.                                 |
| maintenance/<br>assembly                                                                                                        | Lubrication not admitted.                                                                   |
| <b>packaging</b>                                                                                                                | Paper/film bag (single or twice).                                                           |
| <b>sterilization</b>                                                                                                            | Closed and with hooked spring. Not lubricated.                                              |
| recommended<br>classification according<br>to KRINKO/RKI/BfArM<br>guidance (only Germa-<br>ny, with respect to<br>intended use) | Critical B                                                                                  |

|                                                 |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      |
|-------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>Art. no</b>                                  | 8000200                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                              |
| <b>article specification</b>                    | Cartilage Punch Set (see figure 4 and 5)                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                             |
| <b>rinsing volume</b>                           | 50-100 ml (single-use syringe)                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                       |
| <b>brush</b>                                    | inside punch: pipe cleaner 3 mm inside contour punch: pipe cleaner 4 and 6 mm. Outside: standard brush/ tooth brush                                                                                                                                                                                                                                                                                                                                                                                                                                  |
| <b>specific/additional procedure in case of</b> |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      |
| pretreatment                                    | Disassemble into contour punch, inner hole punch, spring, and push-out needle<br>1. Withdraw the push-out needle from the inner hole punch.<br>2. Separate inner hole punch and contour punch. Brush inside and outside.<br>Flush inside.<br>Check for tissue remaining in lumen or on tip (otherwise repeat pre-cleaning or dispose).<br>Plastic Base POM: Brush.<br>Check for damages (e.g. synthetic fibres/particles released from the plastic base) and dispose if damaged.<br>Do not apply specific instrument tray for cleaning/disinfection! |
| manual cleaning/ disinfection                   | Disassembled.<br>Contour punch, inner hole punch, spring, and push-out needle:<br>Brush inside and outside.<br>Flush inside. Plastic Base POM: Brush.<br>Do not apply specific instrument tray for cleaning/disinfection!                                                                                                                                                                                                                                                                                                                            |
| automated cleaning/ disinfection                | Disassembled.<br>Contour punch, inner hole punch, spring, and push-out needle:<br>small pieces basket. Plastic Base POM: Standard tray.<br>Do not apply specific instrument tray for cleaning/disinfection!                                                                                                                                                                                                                                                                                                                                          |
| maintenance/ assembly                           | No assembly prior to sterilization. Lubrication not admitted!                                                                                                                                                                                                                                                                                                                                                                                                                                                                                        |
| <b>packaging</b>                                | Disassembled in specific instrument tray (spring in the hole of the front holder of the Plastic Base POM).<br>Specific instrument tray in paper/film bag (single or twice).                                                                                                                                                                                                                                                                                                                                                                          |

|                                                                                                                |                                                                                                                                                                                                                                                                                      |
|----------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>sterilization</b>                                                                                           | Disassembled in specific instrument tray. Not lubricated.<br>After sterilization/directly prior to application:<br>1. Fit contour punch in place before placing the inner hole punch with spring<br>2. Insert the push-out needle into the inner punch.<br>Lubrication not admitted! |
| recommended classification according to KRINKO/RKI/BfArM guidance (only Germany, with respect to intended use) | Critical B                                                                                                                                                                                                                                                                           |

|                                                                                                                                 |                                   |
|---------------------------------------------------------------------------------------------------------------------------------|-----------------------------------|
| <b>Art. no</b>                                                                                                                  | 8000249 to 8000254                |
| <b>article specification</b>                                                                                                    | Sizer Breathe Implant àWengen     |
| <b>rinsing volume</b>                                                                                                           | -                                 |
| <b>brush</b>                                                                                                                    | standard brush/tooth brush        |
| <b>specific/additional procedure in case of</b>                                                                                 |                                   |
| pretreatment                                                                                                                    | Brush.                            |
| manual cleaning/<br>disinfection                                                                                                | Brush.                            |
| automated<br>cleaning/<br>disinfection                                                                                          | Small pieces basket.              |
| maintenance/<br>assembly                                                                                                        | Lubrication not admitted.         |
| <b>packaging</b>                                                                                                                | Paper/film bag (single or twice). |
| <b>sterilization</b>                                                                                                            | Not lubricated.                   |
| recommended<br>classification according<br>to KRINKO/RKI/BfArM<br>guidance (only Germa-<br>ny, with respect to<br>intended use) | Critical B                        |

|                                                                                                                                 |                                   |
|---------------------------------------------------------------------------------------------------------------------------------|-----------------------------------|
| <b>Art. no</b>                                                                                                                  | 8000555                           |
| <b>article specification</b>                                                                                                    | Sizer Omega Connector             |
| <b>rinsing volume</b>                                                                                                           | -                                 |
| <b>brush</b>                                                                                                                    | standard brush/tooth brush        |
| <b>specific/additional procedure in case of</b>                                                                                 |                                   |
| pretreatment                                                                                                                    | Brush.                            |
| manual cleaning/<br>disinfection                                                                                                | Brush.                            |
| automated<br>cleaning/<br>disinfection                                                                                          | Small pieces basket.              |
| maintenance/<br>assembly                                                                                                        | Lubrication not admitted.         |
| <b>packaging</b>                                                                                                                | Paper/film bag (single or twice). |
| <b>sterilization</b>                                                                                                            | Not lubricated.                   |
| recommended<br>classification according<br>to KRINKO/RKI/BfArM<br>guidance (only Germa-<br>ny, with respect to<br>intended use) | Critical B                        |

## ANNEX TO SPECIFIC ASPECTS

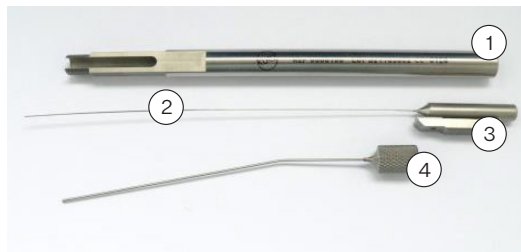
Figure 1:  
**KURZ Meter (8000100)  
/ KURZ Meter Set  
(8000106)**

1. Holder
2. Probe
3. Handle with pusher
4. Tube

Figure 2:  
**Vario-Holder  
(8000135)**

1. Holder, lower part
2. Holder, upper part
3. Pin

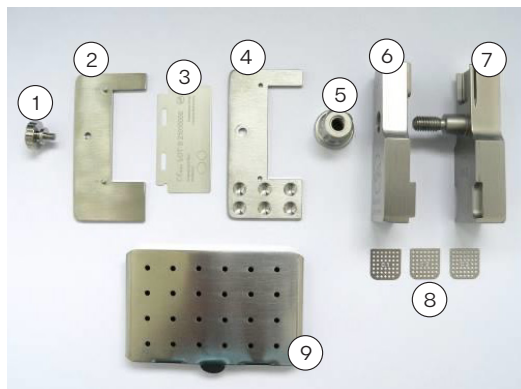
Figure 3:  
**KURZ Precise Cartilage  
Knife Set (8000155)**

1. Screw
2. Blade holder, upper part
3. Single-use blade
4. Blade holder, lower part
5. Nut
6. Cutting block, upper part
7. Cutting block, lower part
8. Distance plates
9. Small Box for accommodating the distance plates

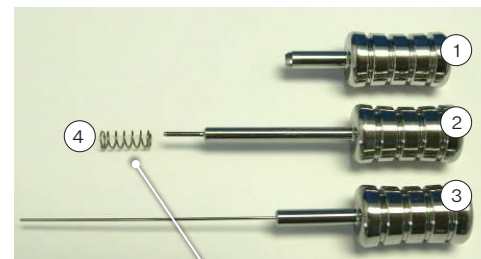
Figure 4:  
**Cartilage Punch Set  
(8000200), Position of  
spring for mounting**

1. Contour punch
2. Inner hole punch
3. Push out punch
4. Spring

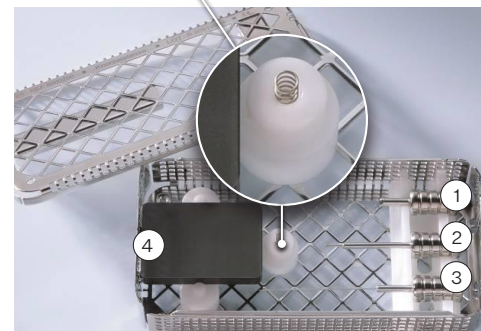
Figure 5:  
**Parts of Cartilage Punch  
Set (8000200) in  
sterilization tray**

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